




PLUS
LABORATORY

PRODUCT SELECTION GUIDE




omni

the most versatile orthotic design available

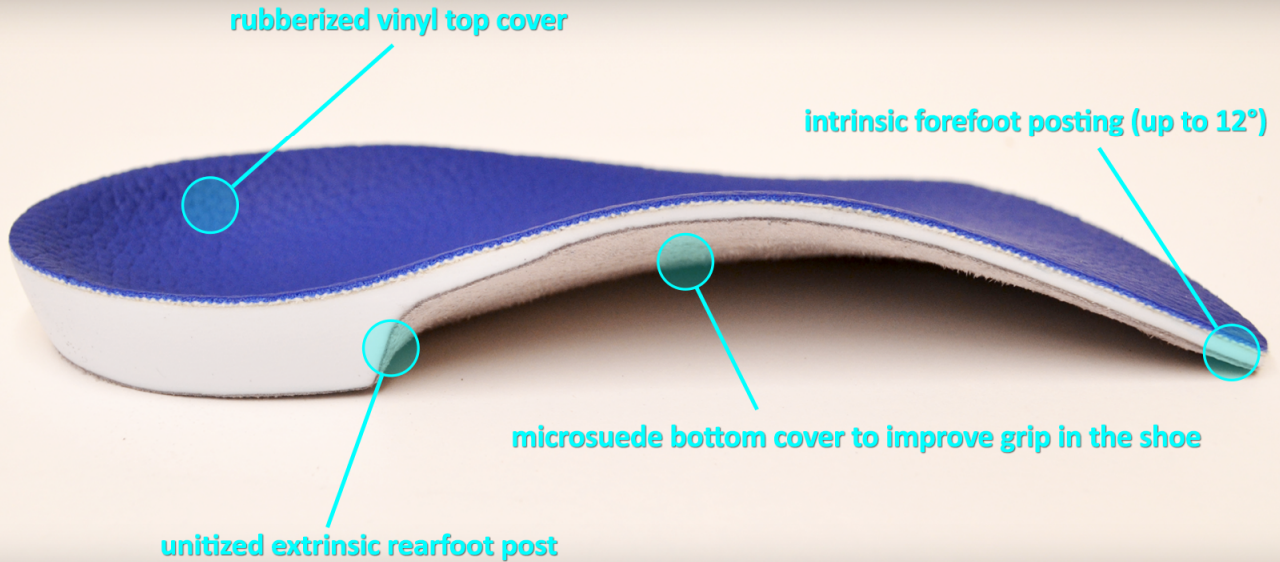

The Omni™ is our most popular design, thanks to a very low bulk profile that is superior to carbon-fiber and all other pressed orthotics. Through its tapered distal edge and thin heel center, both the forefoot and rearfoot achieve minimal bulk while retaining full function and uncompromised biomechanical control.

Customize each Omni by adjusting the Flexibility, Bulk, Posting, Cushioning, the Top Covering, Accommodations, Pads, and Cast Modifications (like Plaster Fill and Cast Inversion) - to achieve your specific goals for that device.

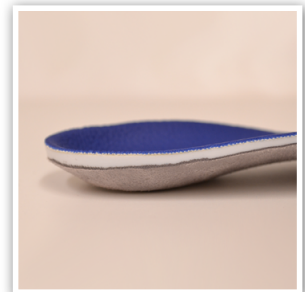

Fashionette

the solution for patient compliance in modern shoes

At Plus Laboratory we understand that the trends in fashion footwear demand a virtual “No Bulk” design. The Fashionette combines a zero-bulk heel, a tapered forefoot edge, and thinned lateral column to create a functional device that can be worn in even the most limiting shoe profiles.

We **guarantee** a perfect shoe fit on all Fashionettes if a shoe or a template of the footbed is provided for our technicians to work from.



Have your patient further tailor their Fashionettes to their shoes and aesthetic preferences by choosing from our broad selection Top Covering colors and fabrics.

VINYL TOP COVERINGS



MICROSUEDES




SOFT LEATHERS



PedFlex

designed to support and protect your growing patients

Our UCBL-type design, the PedFlex provides a deep heel, unitized rearfoot post and medial flange - while maintaining lower bulk in both the rearfoot and forefoot for improved compliance in your pediatric patient.

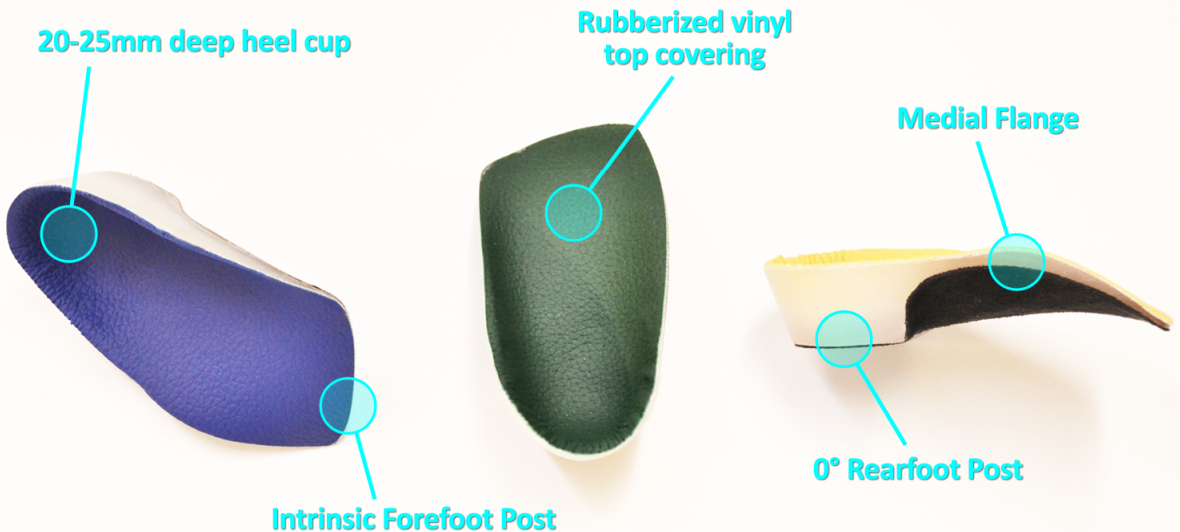

The replacement of outgrown orthotics can be a difficulty for the parents of your young patients. This is why Plus Laboratory was the first orthotics lab to offer an outgrowth protection plan on all our orthotics, purchasable by the parent, that will allow for a single use no-charge fabrication of a new pair of orthotics digitally grown by our CAD/CAM technology to match their child's increase in shoe size.






















Of course, the fabrication of the new orthotics will only be done with your approval, and the orthotics are shipped to your office for dispensing.

DURACCOM

an accommodative device with industry-leading durability


Our modern accommodative design upgrades the traditional leather by utilizing a resilient shell between the soft materials to provide your patients many years of relief without losing its shape.

Complete customization - far beyond what is achievable with Medicare insoles - will be at your fingertips when you choose the DurAccom for your diabetic, elderly, and lesion-prone patients.


| | | | |
|--------------------|---|---|---|
| Durability |  |  |  |
| Low Bulk |  |  |  |
| Fully Customizable |  |  |  |
| Material Options |  |  |  |
| Custom-Molded |  |  |  |
| Lightweight |  |  |  |
| Soft Material |  |  |  |
| | Medicare Insoles | DURACCOM | Traditional Leather |

SUBSTRATE OPTIONS


PORON

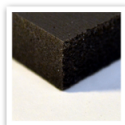


SPONGE

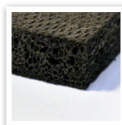

PLASTAZOTE

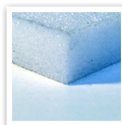

KOREX

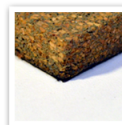

CORK



SPORT PLUS

custom designs for the dynamic lifestyle

ActiveFit



This general purpose, semi-rigid device with rubber rearfoot posts is best suited to the moderately active patient for whom a sport-specific design is unnecessary.



Built for the increased pronation and shock absorption needs of running, with optimum control during the propulsive phase, all while maintaining a low bulk profile to fit modern running shoes.

RunFit



Long Distance



Long Distance Option: Increased shock absorption at the forefoot and heel strike.

FieldFit



The field orthotic is crafted for activities with aggressive lateral acceleration and deceleration. This design also incorporates a narrower profile, to better fit the cleats used for soccer, baseball, football, rugby and lacrosse.

SPORT PLUS

custom designs for the dynamic lifestyle

CourtFit



High Impact



The CourtFit utilizes rubber dynamic rearfoot posts, adjusted lateral controls, and forefoot cushioning to serve the needs of athletes that experience increased forefoot impact and significant lateral forces.

High Impact Option: Full Spenco™ cover to toes, for maximum shock absorption



Your cyclist will appreciate the extra low bulk contouring, slim profile, and custom shoe fitting that comes standard with the CycleFit. This is a rigid device with intrinsic rearfoot posts and inclined Korex extensions to the toes.

CycleFit



GolfFit



The GolfFit is a high flex design, with intrinsic forefoot and intrinsic rearfoot posting for lower bulk, with soft inclined forefoot extensions to the toes.

Orthotic Elements

Unmatched selection. Uncompromised quality.

Top Coverings

Soft Leather
low bulk
durable

Caramel
Blue
Black
Wheat
Espresso
White

Vinyl
rubberized for grip
nylon reinforced

Red
Golden
Green
Saddle
Navy
White
Beige
Gray
Burgundy

Microsuede
lowest bulk
warm

Silver
Camel
Black

Cushioning

Poron
very soft (duro 15)
good shape memory

Plastazote
soft (duro 18)
compresses over time

Spenco
medium (duro 22)
non-slip blue nylon top
good shape memory

EVA
firm (duro 35)
no top covering

Substrates

Poron
very soft (duro 15)

Sponge
soft (duro 17)

Plastazote
firm (duro 30)

Cork
very firm (duro 55)

Korex
very firm (duro 50)
durable

Plastic

White

Natural







TRUE CUSTOMIZATION

from products to service: a perfect fit for your practice

At Plus Laboratory we insure you get the most from your orthotics laboratory, from the beginning of the ordering process to the final product you dispense to your patients.

We specialize in providing the personalized support and customized preferences - in ordering, shipping, and billing - that best serve your practice's individual needs.

Special Requests? That's Our Specialty

-  Shipping to patient or other locations, at no additional cost
-  1-day fabrication rush
-  Email notifications for tracking of high-priority orders
-  Additional fabrication materials (pads, cushioning, etc) available on request for in-office modifications

Build a Design Specific to Your Prescription Style

We will work with you to build and test specialized orthotic designs of your own creation.

In addition, any customized standard orthotic design that fits your needs can become a standard Design option, able to be selected with a single check on all your Plus Laboratory Rx forms.